


# syngo FUSION

1. Aligning and superimposing multi-modality images
2. Open Patient Browser; select the first modality series and load into MPR within the 3D Task-card.
3. Re-open the Patient Browser, select the second modality series.
4. Left click on the Fusion icon; this will add the second modality into the 3D Task-card.
5. Go to the Image sub-task card, and left click on the Registration icon.
6. The Fusion Registration page will appear, left click on the Auto Icon and choose "Precise Registration". Then register, and okay. Choose save registration when this pop-up appears.
7. Left click the Side x Side tool to evaluate the 2 modality series in parallel.

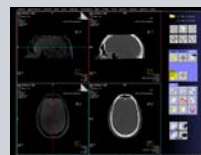

Aligning and superimposing multi-modality images