

Siemens Healthcare Diagnostics, the leading clinical diagnostics company, is committed to providing clinicians with the vital information they need for the accurate diagnosis, treatment, and monitoring of patients. Our comprehensive portfolio of performance-driven systems, unmatched menu offering and IT solutions, in conjunction with highly responsive service, is designed to streamline workflow, enhance operational efficiency, and support improved patient care.

CLINITEK Status, CLINITEK, RAPIDComm, Multistix and all associated marks are trademarks of Siemens Healthcare Diagnostics Inc. All other trademarks and brands are the property of their respective owners.

Product availability may vary from country to country and is subject to varying regulatory requirements. Please contact your local representative for availability.

Siemens Global Headquarters

Siemens AG
Wittelsbacherplatz 2
80333 Muenchen
Germany

Global Siemens Healthcare Headquarters

Siemens AG
Healthcare Sector
Henkestrasse 127
91052 Erlangen
Germany
Telephone: +49 9131 84 - 0
www.siemens.com/healthcare

Global Division

Siemens Healthcare Diagnostics Inc.
1717 Deerfield Road
Deerfield, IL 60015-0778
USA
www.siemens.com/diagnostics

www.siemens.com/diagnostics



How can I effectively manage
point-of-care urinalysis testing?

Flexible connectivity solutions that deliver
testing oversight to point-of-care urinalysis

CLINITEK Status® Connect System

Answers for life.

SIEMENS

Introducing the CLINITEK Status® Connect System,

your solution for better management

CLINITEK Status Connect System brings connectivity and testing oversight to point-of-care urinalysis for improved risk management

Flexible connectivity facilitates integration of urinalysis testing for improved documentation and billing

Connect to your hospital's LIS/HIS or point-of-care data management software

Store data from multiple testing locations in a central database

Automatically transmit patient data and QC results

Flexible connectivity via serial, ethernet or wireless connection and supports HL7 and POCT1-A2 communication protocol



Advanced operator management provides an additional level of security to improve risk management

Stores up to 700 operator IDs per analyzer

Prevents unauthorized use based on customizable settings

Programmable to grant specific levels of access for each operator



Expanded Quality Control helps managers and supervisors drive compliance across testing sites

Define specific QC protocols for operators to follow

QC lock-out if protocol is not followed or QC test fails

New security features and quality checks improve data accuracy

Checks each strip automatically for humidity to prevent reporting of falsely elevated results (available on strips with leukocyte pad)

Detects common sample interferences in urine and provides printed notation for improved results interpretation* (availability dependent on strip type)



Comprehensive urine test menu – including hCG pregnancy testing

Enter data using optional barcode reader to eliminate manual steps and reduce transcription errors

*Feature not available in the U.S.

Want testing oversight and control of multiple point-of-care urinalysis analyzers and operators?

CLINITEK Status Connect System improves risk management when interfaced with point-of-care data management solutions

Manage your point-of-care urinalysis testing program more efficiently

- Gain centralized control and oversight for point-of-care urinalysis analyzers in multiple locations
- Maintain up-to-date certification records of operators and remotely restrict operator access if needed
- Transmit patient results electronically, eliminating the need for manual documentation
- Generate QC documentation and reports whenever needed

CLINITEK Status Connect System interfaces to leading data management solutions:

- Siemens' dedicated RAPIDComm® Data Management and Connectivity Solution
- Conworx Technology POCcelerator™*
- MAS RALS®-Plus*
- TELCOR Quick-Linc®*
- For others, contact your data management vendor




*Interfaces in development.



Receive urinalysis results from offsite physician's office.

Connecting to data managers provides visibility to urinalysis testing in decentralized sites

Bi-directional communication only when using a data management middleware solution

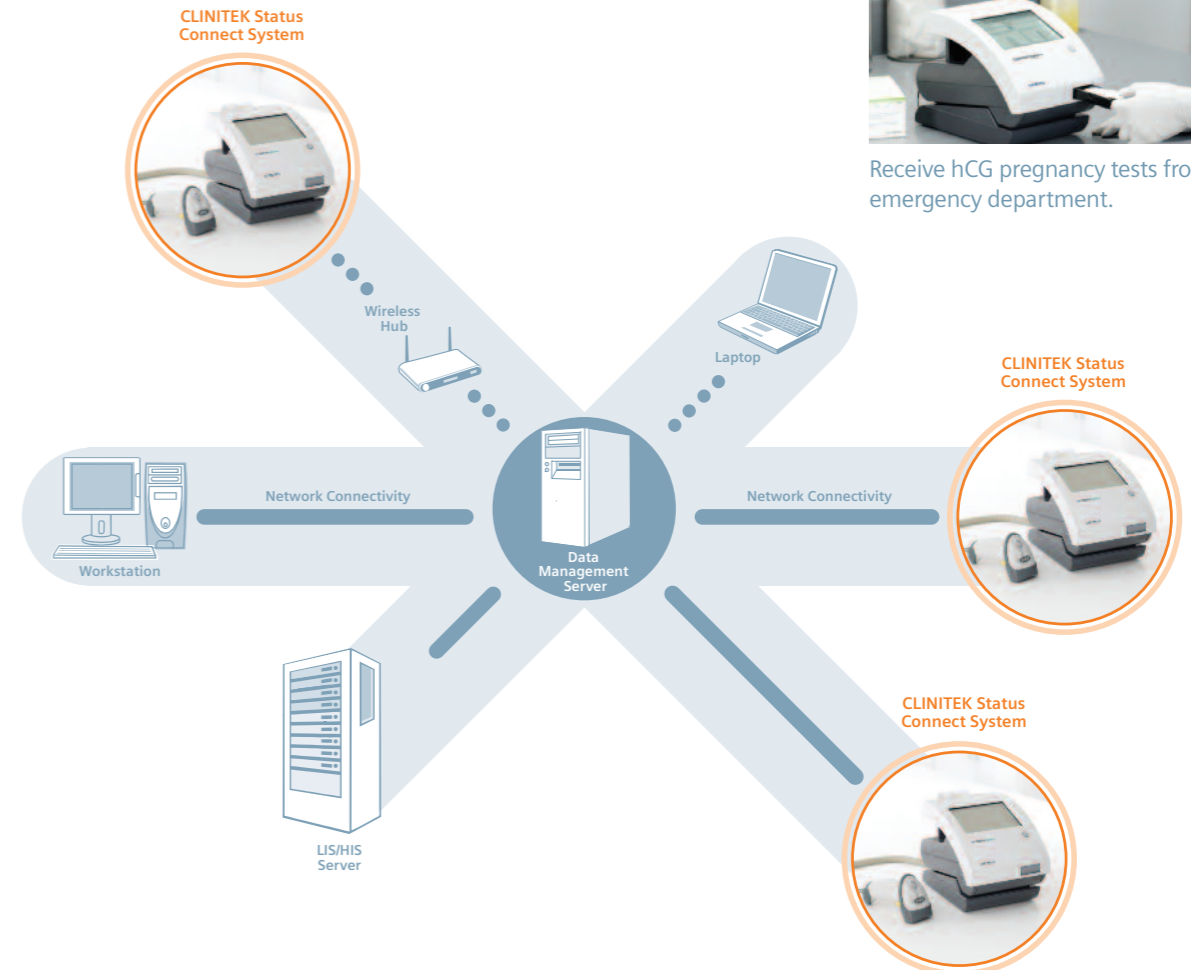

Receive hCG pregnancy tests from emergency department.

RAPIDComm Data Management and Connectivity Solution – a single-vendor option from Siemens



RAPIDComm Data Management and Connectivity Solution
Version 3.0 – Remotely manage multiple CLINITEK Status
Connect Systems from a central location

A complete, efficient urinalysis solution from one source:
Why RAPIDComm Data Management and Connectivity Solution?

Manage data from multiple analyzers in one central location

- Review QC data and generate customized reports
- Query hospital system for patient demographics
- Manage operator data and download operator lists to analyzers

Remotely monitor and control multiple analyzers in multiple locations without physically visiting sites

- Monitor analyzer status and events
- QC lockout – prevent testing if QC regimen is not followed or QC fails
- Perform operator lockout if recertification is overdue
- Automatically recertify operators
- Customizable set-up configuration can be downloaded to all analyzers

CLINITEK – a legacy of innovation in urinalysis
testing from Siemens Healthcare Diagnostics

Trust and confidence based on a heritage of delivering
reliable solutions

- Most widely used system worldwide; over 60,000 CLINITEK analyzers shipped
- CLINITEK and Multistix® – the brands clinicians know and trust
- CLIA-waived and CE-marked

Siemens leads the way in urinalysis
testing with world-class service

To learn more about this exciting new offering,
contact your local Siemens Healthcare Representative
or visit our web site at www.siemens.com/urinalysis
or www.usa.siemens.com/urinalysis.