

Abdominal Imaging

REAL SOLUTIONS

Understanding challenges
Providing answers
Improving outcomes



The ACUSON Antares™ ultrasound system, premium edition is ideally suited for abdominal imaging. Unique features such as Hanafy lens transducers and Clarify VE technology provide superior image quality on even the most technically difficult-to-image patients. TEQ technology improves image quality while enhancing clinical workflow. With its compact architecture and advanced ErgoDynamic™ imaging system design, the Antares system is the ideal clinical solution for today's abdominal imaging challenges.

Highlights

Proven technically difficult-to-image patient performance:

- 30 cm display depth
- Extend imaging technology
- Tissue Harmonic Imaging
- Hanafy Lens transducer technology

Advanced SieClear™ spatial compounding*:

This real-time compounding technique applies industry-first 13 lines of sight at greater steering angles to improve contrast resolution and border detection. Siemens' Dynamic TCE algorithm simultaneously achieves speckle reduction, contrast enhancement and improvement in coherence of anatomical structures. Tissue Stabilization reduces the artifacts seen in compounded images, yielding unprecedented image clarity.

Clarify™ vascular enhancement technology*:

Reduces noise and artifacts to offer exceptional contrast resolution and clear delineation of micro- and macrovasculature in abdominal organs and vessels.

Cadence™ contrast pulse sequencing technology* **:

Enhances sensitivity and specificity to contrast-agent signals for excellent lesion detection and characterization.

TEQ™ ultrasound technology:

Provides single-button optimization for 2D and spectral Doppler imaging, optimizing images for consistent reproducible image quality and reduced operator strain.

Extend imaging technology*:


Improves penetration at depth in 2D, Color and PW Doppler to enhance image quality in technically difficult-to-image patients.

ACUSON Antares ultrasound system, premium edition
The ideal platform for today's abdominal imaging challenges




Advanced SieClear spatial compounding with Dynamic TCE technology

- Provides excellent speckle reduction, contrast enhancement and improved coherence of anatomical structures as seen in this image of a liver and kidney.

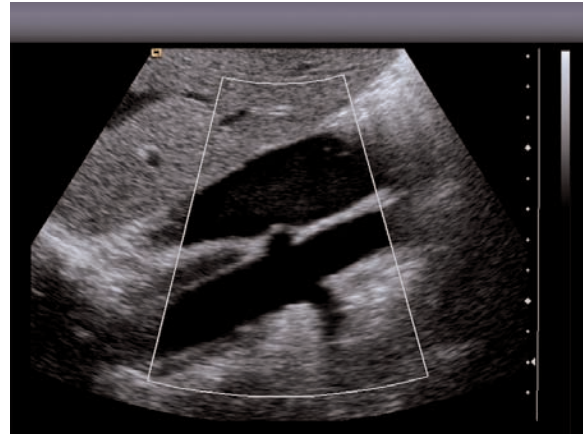


Clarify vascular enhancement technology

- Removes noise and artifacts to clearly demonstrate the aorta with renal arteries and inferior vena cava.

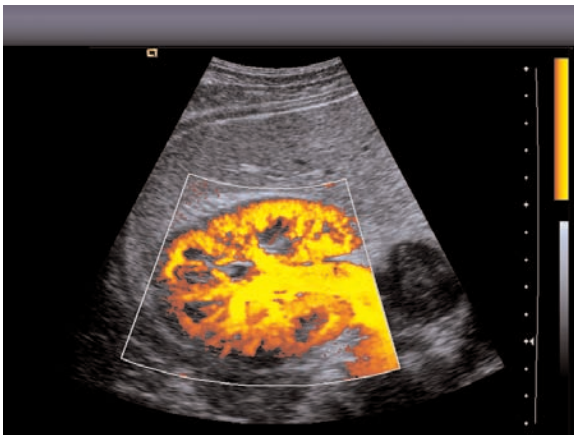

Power Doppler imaging

- Demonstrates excellent sensitivity and spatial resolution as seen in this image of normal renal vasculature.

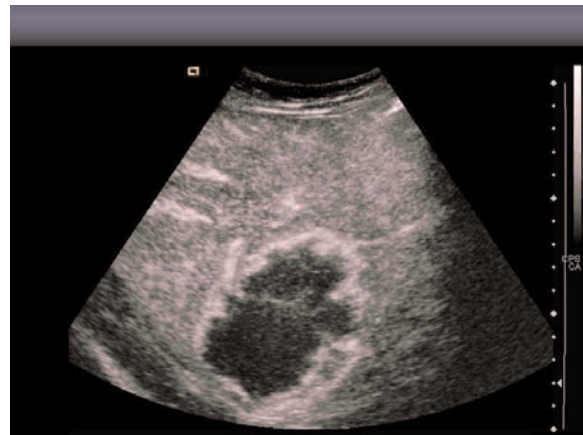

Cadence contrast pulse sequencing technology**

- Delivers unparalleled sensitivity and specificity to contrast imaging agents as demonstrated in this hemangioma of the liver.

*Available on ACUSON Antares, premium edition only.

**At the time of publication, the US Food and Drug Administration cleared ultrasound contrast agents only for use in LVO. Check the current regulation for the country in which you are using the system for contrast agent clearance.

© 02.2007, Siemens AG
 Order No.
 A91US-47-1C-4A00
 Printed in USA
 WS 0207 5.0

Headquarters
 Siemens Medical Solutions USA, Inc.
 51 Valley Stream Parkway
 Malvern, PA 19355-1406 USA
 Telephone: +1-888-826-9702
www.usa.siemens.com/medical

Europe: +49 9131 84-0
 Asia Pacific: +65 6341 0990
 Latin America: +1-786-845-0697

USA
 Siemens Medical Solutions USA, Inc.
 Ultrasound Division
 1230 Shorebird Way
 P.O. Box 7393
 Mountain View, CA 94039-7393 USA
 Telephone: +1-888-826-9702